

Klassik-Konso Knee Joint



Features

Sizes available for the **smallest to large** patients

Joint flexes to **140°**

Fabricated from **high-strength aluminum**

+ Adjustable flexion angle from **0° to 40°** in 5° intervals

Also available with a **straight lower bar** for easier separation of the thigh from calf section

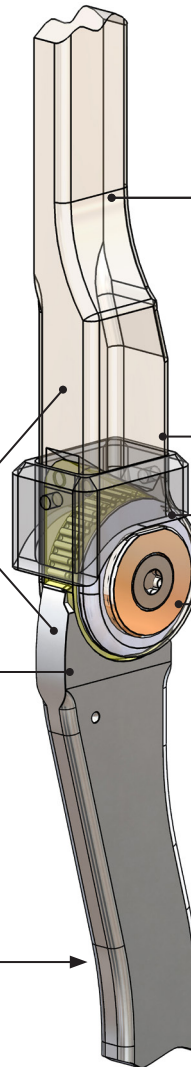


Fig. left

Standard

Fully assembled as one item for clear cost estimates and uncomplicated ordering

Integrated upright bars ensure **optimal transfer of force**

Upper section **CNC-machined from one piece** of material for increased stability

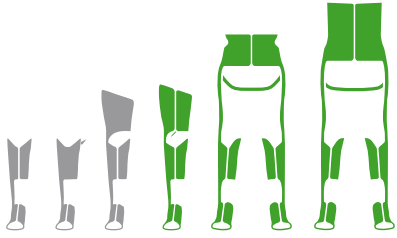
Precise assembly calibrated to 1/100mm with self-lubricating washers for smooth operation free of play

Maintenance-free stainless steel joint assembly with bronze bushing

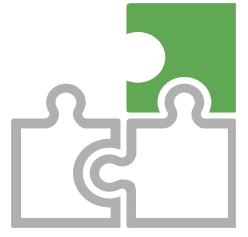
Unrivalled **20-month guarantee**

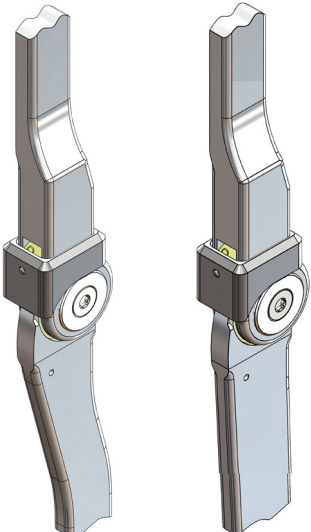
Engraved **serial number** for traceability

Clearly illustrated instructions for use


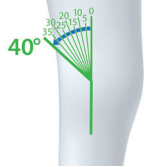
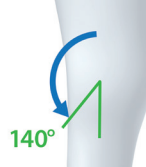

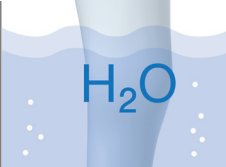


We recommend using our configurator to make your selection.
The configurator considers, among other factors, weight to height ratios.



Weight Classification (body weight)			Joint			
Article No.	KAFO	HKAFO/RGO		Width	Thickness	Drop Lock Thickness
Klassik-Konso-22	up to 10 kg	up to 10 kg		22 mm	10,5 mm	15,0 mm
Klassik-Konso-25	up to 20 kg	up to 15 kg		25 mm	13,0 mm	16,5 mm
Klassik-Konso-28	up to 30 kg	up to 20 kg		28 mm	15,0 mm	19,5 mm
Klassik-Konso-30	up to 35 kg	up to 30 kg		30 mm	16,0 mm	20,5 mm
Klassik-Konso-32	up to 50 kg	up to 40 kg		32 mm	17,0 mm	21,0 mm
Klassik-Konso-36	up to 70 kg	up to 50 kg		34 mm	17,0 mm	21,0 mm
Klassik-Konso-40	up to 85 kg	up to 65 kg		40 mm	21,0 mm	26,0 mm
Klassik-Konso-40-Plus	up to 85 kg	up to 85 kg	40 mm	22,0 mm	28,0 mm	

Klassik-Konso contoured | Klassik-Konso straight

								
Locked	Extension	Flexion	Posterior Offset	Modular	Force Options	No Grinding Necessary	Noise Dampening	Water Resistant